

# "V" MIXER

FTLMV




The laboratory "V" mixers are used to homogenize and mix dusty products. They have two inlets with silicone lids and one outlet for the resulting product. The optimal capacity to homogenize is 50% of the total volume mixer.

The body is made of st. steel AISI 304 with a polished mirror in both internal and external sides. It also can be manufactured in AISI 316.

It works with a single-phase motor and has a timer to control the mix time.

This "V" container is made without corners and it is brilliant polished (inside and outside) to avoid dirty and make an easy cleaning.

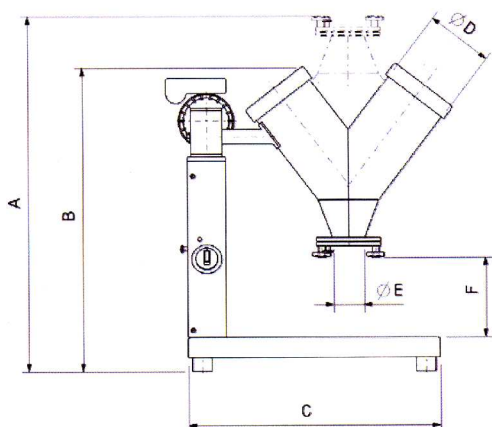

(mm)	2 l.	4 l.	8 l.
A	500	705	805
B	445	605	700
C	415	500	600
D	104	133	133
E	60	60	60
F	128	158	200

## Specifications and ordering information

- laboratory "V" mixer 2 litres
- laboratory "V" mixer 4 litres
- laboratory "V" mixer 8 litres

.....  
 .....  
 .....

Cat. No. FTLMV-0002  
 Cat. No. FTLMV-0004  
 Cat. No. FTLMV-0008

### Technical specifications:

Power: 230Volt - 50/60Hz - 0.12CV (0.25CV for FTLMV-0008)



www.socachim.com

Av. du Capricorne 149, 1200 Brussels/Belgium

TEL: +32 (0)2 762 77 12 FAX: +32 (0)2 762 55 07

socachim@socachim.com